



Product Catalog



3925B


Name	Number	Weight		Description
BASE SLED	EN-BU-002	lb 38	kg 17.2	Carbon-fiber Base Sled serves as a sturdy foundation platform where radiation protection components attach. This sled comes equipped with a built-in CPR board.
MATTRESS TOPPER	EN-MT-001	8.5	3.9	For maximum patient comfort, the 1.5" foam mattress topper is secured on top of the Base Sled.
RIGHT & LEFT ARM BOARDS WITH INTEGRATED FLIP SHIELD	EN-AB-R01 EN-AB-L01	4.5	2.0	Carbon-fiber Arm Board attaches to the Base Sled to offer reliable support for radial procedures. The Arm Board is equipped with innovated flip shields, effectively blocking scatter radiation.
LEFT & RIGHT HIP FLIP SHIELDS	EN-HF-R01 EN-HF-L01	2	0.9	Carbon-fiber Flip Shields attach to the right and left side of the Base Sled rail to provide radiation protection above the table at the hip of the patient.
CLEAR SPOT SHIELD	EN-CS-001	12	5.4	Transparent Acrylic Shield with adjustable arm post provides manual positional protection from the C-arm above the table.

Name	Number	Weight		Description
RIGHT & LEFT SIDE HANGING FLEX SHIELDS	EN-SS-R01 EN-SS-L01	lb 14	kg 6.4	Hanging Shielding System is supported by internal stays, allowing smooth movement of the C-arm without compromising safety. Seamlessly attach to the Base Sled rail and protects from scatter radiation below the table on the sides.
HANGING HEAD FLEX SHIELD	EN-HS-001	8.5	3.9	Hanging Shielding System protects from scatter radiation below the table at the head. It is securely attached to the rail that wraps around the Base Sled head.
STRUCTURAL HEAD SHIELD	EN-HB-001	13.5	6.1	Fits on top of the mattress to provide increased protection from scatter radiation above the table at the head.
WORKBENCH	6" 12" 36"	Varies by Model	Varies by Model	Incorporates flip shielding to block scatter radiation from the lower extremity. Manufactured in various configurations to accommodate different physician-technician workflows.

Name	Number	Weight		Description
MATTRESS EXTENDER	EN-ME-001	lb R: 4 L: 5.5	kg R: 1.8 L: 2.5	This mattress is designed to extend from the end of the Base Sled and Mattress Topper to the end of the operating table, creating a seamless and even work surface. Comes in regular or long sizes.
EGG PAD	EN-EP-001	3	1.4	Reusable radiation pad placed on patient to protect from patient scatter. Effective replacement of RadPad.
EGG PAD STERILE DRAPE	EN-PD-001	0		Sterile Drape that is designed to fit securely and easily over the Egg Pad.
FLEXIBLE ANESTHESIA HOOP	EN-AH-001	2.25	1.0	This anesthesia hoop shapes and molds easily over the patient and table to allow for optimal workflow and drape positioning.
CLEAR SPOT SHIELD STERILE DRAPE	EN-SD-001	0		Transparent Sterile Drape that is designed to fit securely and easily over the Clear Spot Shield.

EGGNEST™


XR

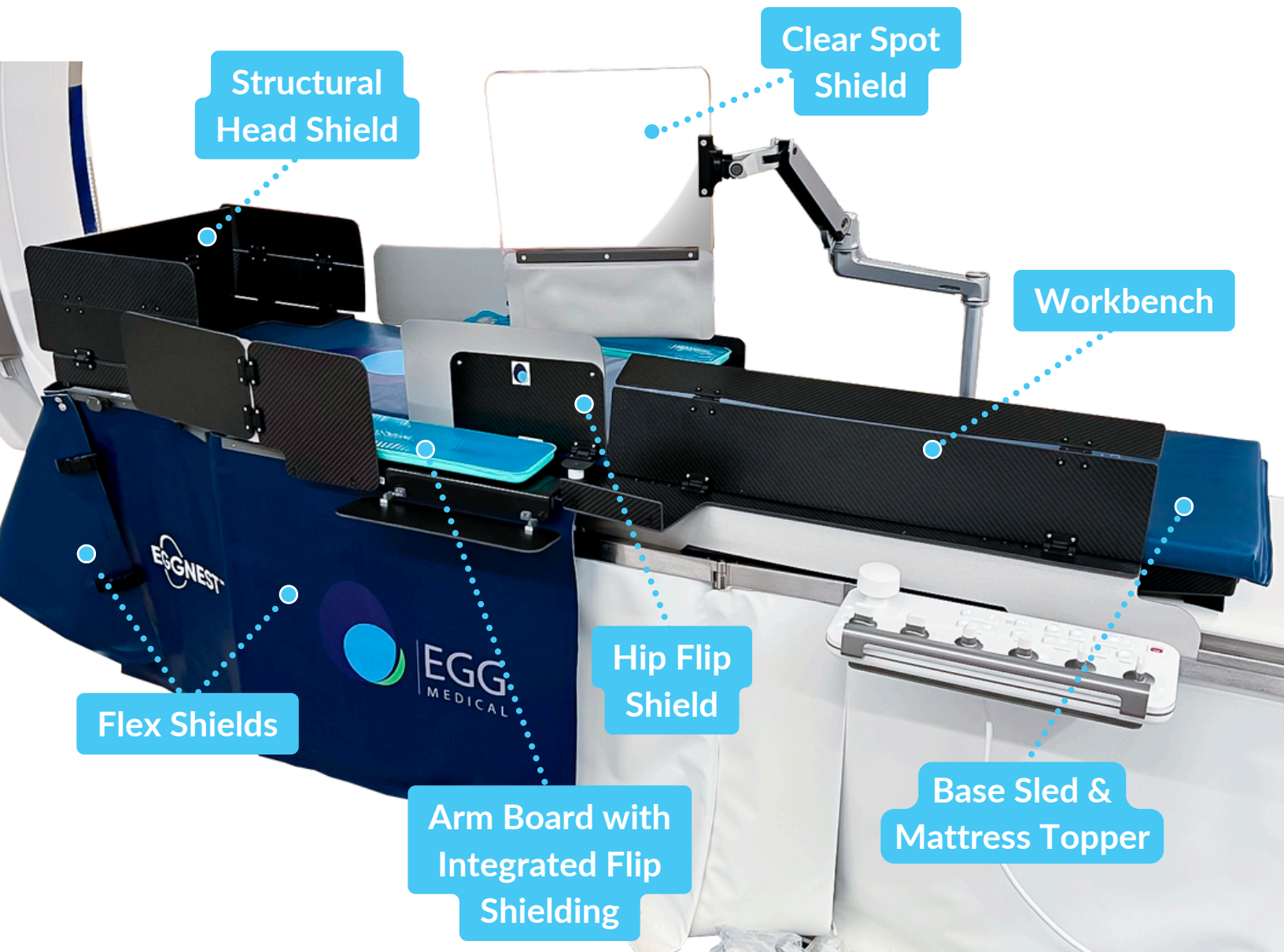

Order Your Components Today:

 eggmedical.com

 info@eggmedical.com

 (612) 916-6616



3925B